

仅限科研

版本号：202301

β-Actin Rabbit Polyclonal antibody

产品内容


名称	货号	规格	保存
β-Actin Rabbit Polyclonal antibody	SDP0043	50ul	-20℃

产品描述

产品名称	β-Actin Rabbit Polyclonal antibody
宿主	Rabbit
靶标种属	Human、Mouse、Rat
应用	WB: 1:500-5000; ICC/IF: 1:10; ELISA: 1:500-2000; IHC: 1:100-1:200
产品形式	Liquid
存储溶液	0.01M TBS(pH7.4) with 1% BSA, 0.03% Proclin300 and 50%Glycerol
浓度	1 mg/mL
纯度	Immunogen affinity purified
纯化说明	This antibody was isolated by affinity chromatography using antigen coupled to agarose beads and conjugated to HRP.
克隆	Polyclonal
同种型	IgG
抗原	Gene ID: 60; Swiss Prot: P60709; Synthetic MAP peptide derived from human beta-Actin:1-200/375
别名	Beta Actin; beta-Actin; ACTB; Actin cytoplasmic 1; Actin, beta; Beta actin; betacytoskeletal actin; A X actin like protein; ACTB; Actin cytoplasmic 1; alpha sarcomeric Actin; Actx;Beta cytoskeletal actin; Melanoma X actin; PS1TP5BP1; ACTB_HUMAN. β actin; β actin
大小	42kd




产品应用

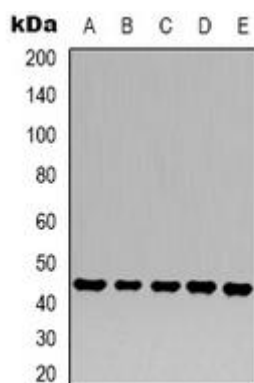

Western blot analysis of Beta-actin expression in 293T (A), HeLa (B), mouse kidney (C), rat heart (D), rat brain (E) whole cell lysates

Primary: β -Actin Rabbit Polyclonal antibody (Cat:SDP0043) at 1/1000 dilution;

Secondary:HRP Goat Anti-Rabbit IgG (H+L) Antibody (cat:SDP0030) at 1/2000 dilution;

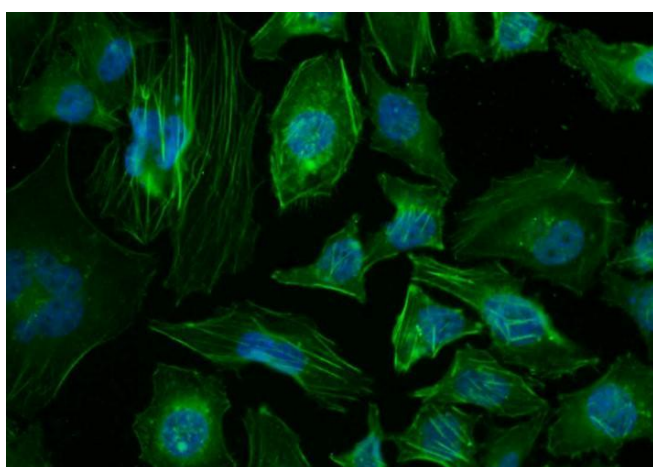

Tissue/cell: HeLa cell;

Primary: β -Actin Rabbit Polyclonal antibody (Cat:SDP0043)

at 1/10 dilution;

Secondary:HyperFluor 488 Goat Anti-Rabbit IgG (H+L) Antibody (cat:SDI0010) ;

DAPI(cat: SDC0046) was used to stain the cell nuclei

